Deep-Tissue

Percussion Therapeutic Massager

The Brand Used By Professionals® Since 1919

© 2011 Wahl Clipper Corporation • Part No. 93342-900 • Printed in China
www.wahlmassagers.com

IMPORTANT SAFEGUARDS

When using an electrical appliance, basic precautions should always be followed, including the following:

READ ALL INSTRUCTIONS AND SAFEGUARDS BEFORE USING.

DANGER

To reduce the risk of death or electric shock:

1. Do not reach for an appliance that has fallen into water. Unplug immediately.
2. Do not use while bathing or in the shower.
3. Do not place or store appliance where it can fall or be pulled into a tub or sink. Do not place in or drop into water or other liquid.
4. Always unplug this appliance from the electrical outlet immediately after using and before cleaning.
5. Do not use massager without attachments.


WARNING

To reduce the risk of burns, fire, electric shock or injury to persons:

1. This appliance should never be left unattended when plugged in. Unplug from outlet when not in use, and before putting on or taking off parts.
2. Do not operate this appliance under a blanket or pillow. Excessive heating can occur and cause fire, electric shock or injury to persons.


1

MESSAGE ATTACHMENTS

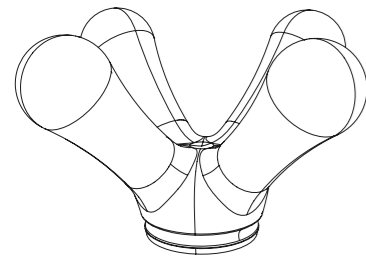
FOUR FINGER FLEX ATTACHMENT

Soft, finger-like nodes replicate a human touch massage.

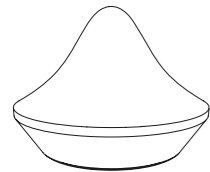
ACCUPOINT ATTACHMENT

Targets deep into tissue. Great for pinpointing concentrated areas.

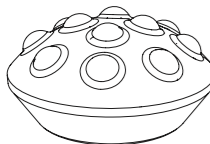RAISED BUMP ATTACHMENT

Wide surface for an all-over muscle massage.

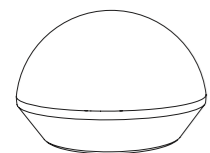DEEP MUSCLE ATTACHMENT

Thumps smoothly, creating a gentle, relaxing massage for all-over general body use.

5

CHANGING ATTACHMENTS

To remove an attachment head, turn the head counter-clockwise. To insert an attachment head, place the head on the stem and turn it clockwise until secure.

CLEANING/STORAGE

Always unplug massager before cleaning.

To clean, simply wipe unit with a dry, soft cloth. Do not immerse in water or use any harsh cleaners or solvents. When storing unit, make sure cord is not twisted or kinked.

6

3. Close supervision is necessary when this appliance is used by, on, or near children, invalids or disabled persons.
4. This massager contains an electro magnetic motor. To avoid risk of injury, it must not be used within six inches of an implantable cardiac device.
5. Use this appliance only for its intended use, as described in this manual. Do not use attachments not recommended by the manufacturer.
6. Never operate this appliance if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Return to Wahl for examination and repair.
7. Do not carry this appliance by supply cord or use cord as a handle.
8. Keep the cord away from heated surfaces.
9. Never drop or insert any object into any opening on this appliance.
10. Do not use outdoors.
11. Do not operate where aerosol (spray) products are being used or where oxygen is being administered.
12. In case of pregnancy, diabetes or illness, consult your physician before using this massager.
13. Do not use massager on:
 - a. Sensitive skin, swollen or inflamed areas, in the presence of poor circulation, where skin eruptions are present, or in the presence of unexplained calf or abdominal pain.

2

MESSAGE TECHNIQUES

SHOULDERS AND NECK —

Massage each shoulder toward back of neck to point just behind and below ear. Then back to shoulder. Repeat several times.

ARMS —

Start with a gentle massage of the palms and fingers. Next, using a circular motion, move upward to forearms. Then up to biceps and shoulders.

LEGS —

If desired, massage bottoms of feet. Then, using a circular motion, massage each leg from calf on up and over thigh.

CHEST —

Massage from center of chest outward, and up toward each shoulder.

BACK —

Massage from top of spine to bottom, and from center outward.

ABDOMEN —

Relax abdominal muscles by massaging from center outward, using a gentle, circular motion.

JOINTS —

Relieve tired, aching knees and elbows with smooth, circular strokes.

FACE —

Gently massage from the nose outward, taking care to avoid contact with the eye.


SCALP —

Start at the back of the neck and massage up to the top of the head. Using this up and down motion, massage across the scalp.

7

- b. Varicose veins.
 - c. The front of the neck.
 - d. Genital areas.
 - e. Sleeping or unconscious persons.
14. When using a massage attachment, be sure to secure the attachment fully onto the massager and keep it on while massaging.
 15. Do not use this massager for more than 20 minutes at one time.
 16. To disconnect, turn unit to the OFF position, then remove plug from outlet.

SAVE THESE INSTRUCTIONS For Household Use Only

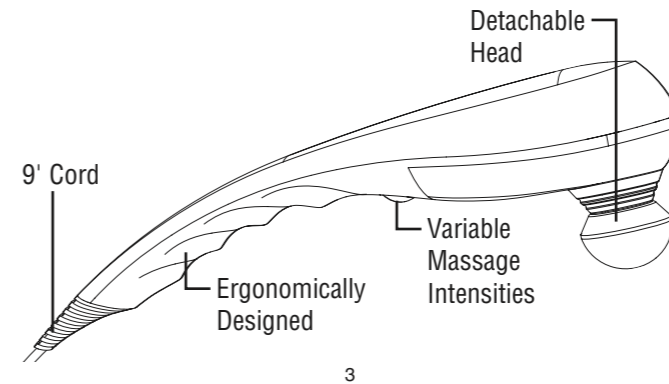3

SERVICE AND REPAIR

This unit has no user-serviceable parts. If you have problems with your massager, **do not return it to store where purchased.** Send it prepaid with a note explaining the problem to:

Wahl Clipper Corporation
Attention: Repair Department
3001 Locust Street
Sterling, Illinois 61081
Telephone: (800) 767-9245

For products purchased outside the United States, contact your Wahl distributor.

For more information or Product Registration, visit www.wahlmassagers.com

8

OPERATING INSTRUCTIONS

To reduce the risk of electric shock, this appliance has a polarized plug (one blade is wider than the other). This plug will fit in a polarized outlet only one way. If the plug does not fit fully into the outlet, reverse the plug. If it still does not fit, contact a qualified electrician to install the proper outlet. Do not change the plug in any way.

1. With the switch in the OFF position, plug the power cord into an AC outlet. (Make sure the voltage of power supply is the same as that indicated on massager.)
2. Place the desired massage attachment on the massager. The attachments easily twist on and off.
3. Variable speed allows you to personalize your massage. For best results, begin at the lowest speed setting, and gradually increase to the desired massaging speed by rolling the dial toward the top of the massager. Roll the dial toward the bottom of the massager to decrease the speed.
4. Apply the massager to the targeted area, using continual movement. A circular motion is often most effective. (See MESSAGE TECHNIQUES.)

4

2-YEAR LIMITED WARRANTY

Valid in U.S.A. Only

If your Wahl product fails to operate satisfactorily within two (2) years from the date of original purchase or receipt as a gift, Wahl will, at our option, repair or replace it and return it to you FREE OF CHARGE. Please mail your product prepaid to Wahl Clipper Corporation, 3001 Locust Street, Sterling, Illinois 61081, Attention: Repair Department, with a note describing the defect and the date of original purchase or receipt as a gift (if in warranty). Products out of warranty will be repaired at our standard repair charges.

This warranty is void and no repair or replacement will be made under this warranty or otherwise if:

1. The product is modified in any manner or repaired by anyone other than Wahl Clipper Corporation or Wahl authorized service center.
2. The product has been used commercially or subjected to unreasonable use.

This warranty gives you specific legal rights and you may also have other rights which may vary from state to state. No other written express warranty is given. Any implied warranty, including any warranty of merchantability which may arise from purchase or use is limited to the years period provided in the express warranty. Some states may not allow such limitation, so it may not apply to you. No responsibility is assumed for incidental or consequential damages of any type or for the use of unauthorized attachments. Some states may not allow the exclusion or limitation of incidental or consequential damages, so the limitation or exclusion may not apply to you.

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WAHL GmbH Roggenbachweg 9 D-78089 Unterkirnach Germany Phone: 49-7721-8060 Fax: 49-7721-806114